



DOCUMENT/RECORD

Identifier

Author

Benjamin Isaacson

Pages:

(including this page) 129

Signed By

Date/Time (GMT)

05/22/2014 07:31:59 PM

05/22/2014 07:33:24 PM

05/22/2014 07:45:56 PM

[illegible]



Medtronic

Nexus-D API Interface Specification

CAUTION – Investigational device. Limited by Federal (or United States) law to investigational use.

1. Scope

This document is intended to be a supplement to the User Manual, Nexus-D API. It contains the lower level technical details that an engineer intending to programmatically interface to the Nexus-D API would be interested in. It does not contain any details about the Nexus-D System Interface.

2. References

Document	Title
NRP1088-36356	User Manual, Nexus-D API
NRP1088-35790	Interface Specification, Nexus-D System

Nexus-D API Documentation

Package Summary		Page
mdt.neuro.nexus		3
mdt.neuro.nexus.commands		54
mdt.neuro.nexus.data		103
mdt.neuro.nexus.support		127

Package mdt.neuro.nexus

Interface Summary		Page
IDisposable	The Interface IDisposable.	19

Class Summary		Page
ApplicationCommand	The Class ApplicationCommand.	4
CommandAccessor	The Class CommandAccessor.	9
DataPacketAccessor	The Class DataPacketAccessor.	17
NexusInstrument	The Class NexusInstrument.	20
ObservableConnection	The Class ObservableConnection.	30
ParameterizedRunnable	The Class ParameterizedRunnable.	34
SerialConnection	The Class Connection.	39
SerialPortErrors	The Class SerialPortErrors.	42
ThreadedNexusInstrument	The Class ThreadedNexusInstrument.	43
ThreadManager	The Class ThreadManager.	50
Util	The Class Util.	52

Enum Summary		Page
ConnectionStatus	The Enum ConnectionStatus.	15
ResponseCodes	The Enum ResponseCodes.	36

Class ApplicationCommand

mdt.neuro.nexus

java.lang.Object

└ mdt.neuro.nexus.ApplicationCommand

Direct Known Subclasses:

CycleStim, DecAmplitude, DecFrequency, DecPulseWidth, GetInsInfo, GetNexusStatus, GetRealtimeData, IncAmplitude, IncFrequency, IncPulseWidth, PulseStim, ResetCycle, RestoreClinicianSettings, SendTrigger, SetActiveGroup, SetNexusConfiguration, StartRealTime, StartSensing, StopRealTime, StopSensing, TherapyOff, TherapyOn

```
public class ApplicationCommand
extends Object
```

The Class ApplicationCommand. This class represents a command to be sent to the Nexus-2 system.

Field Summary		Page
protected List<Byte>	commandBytes The command data (command data plus command code).	6
protected short	commandCode The command code.	6
protected boolean	crcValid The crc valid.	6
protected List<Byte>	request The request.	7
protected List<Byte>	response The response.	6
protected short	responseCode The response code.	6
protected List<Byte>	responseData The response data (not including crc, command code, or response code).	6
protected int	timeout The timeout.	6

Constructor Summary		Page
protected	ApplicationCommand (List<Byte> bytes, int timeoutMs) Instantiates a new application command.	7
protected	ApplicationCommand (short code, int timeoutMs) Instantiates a new application command.	7

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	7
protected List<Byte>	getRequestData () Gets the request data (does not include CRC).	7
protected List<Byte>	getResponseBytes () Gets the response bytes that comprise this payload (including CRC).	8
protected List<Byte>	getResponseData () Gets the response data (does not include CRC).	7

Field Detail

commandCode

protected short **commandCode**

The command code.

commandBytes

protected List<Byte> **commandBytes**

The command data (command data plus command code).

responseCode

protected short **responseCode**

The response code.

responseData

protected List<Byte> **responseData**

The response data (not including crc, command code, or response code).

crcValid

protected boolean **crcValid**

The crc valid.

timeout

protected int **timeout**

The timeout.

response

protected final List<Byte> **response**

The response.

request

protected final `List<Byte>` **request**

The request.

Constructor Detail

ApplicationCommand

```
protected ApplicationCommand(short code,  
                               int timeoutMs)
```

Instantiates a new application command.

Parameters:

code - the command code
timeoutMs - the command timeout (ms)

ApplicationCommand

```
protected ApplicationCommand(List<Byte> bytes,  
                               int timeoutMs)
```

Instantiates a new application command.

Parameters:

bytes - the bytes to be sent
timeoutMs - the timeout ms

Method Detail

getRequestData

```
protected List<Byte> getRequestData()
```

Gets the request data (does not include CRC).

Returns:

the request data

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

getResponseData

```
protected List<Byte> getResponseData()
```

Gets the response data (does not include CRC).

Returns:
the response data

getResponseBytes

protected `List<Byte>` **getResponseBytes**()

Gets the response bytes that comprise this payload (including CRC).

Returns:
the response bytes

Class CommandAccessor

mdt.neuro.nexus

java.lang.Object

└ mdt.neuro.nexus.CommandAccessor

```
abstract public class CommandAccessor
extends Object
```

The Class CommandAccessor. This class allows friend access to the creation of commands by internal API classes not in the commands package.

Constructor Summary	Page
CommandAccessor ()	10

Method Summary	Page
protected abstract getCycleStimCommand () CycleStim Gets the cycle stim command.	10
protected abstract getDecAmplitudeCommand (byte progNum, byte numSteps, byte repeats) DecAmplitude Gets the dec amplitude command.	10
protected abstract getDecFrequencyCommand (byte groupNum) DecFrequency Gets the dec frequency command.	11
protected abstract getDecPulseWidthCommand (byte progNum) DecPulseWidth Gets the dec pulse width command.	11
protected abstract getGetInsInfoCommand () GetInsInfo Gets the gets the ins info command.	11
protected abstract getGetNexusStatusCommand () GetNexusStatus Gets the gets the nexus status command.	11
protected abstract getGetRealTimeDataCommand () GetRealtimeData Gets the gets the real time data command.	11
protected abstract getIncAmplitudeCommand (byte progNum, byte numSteps, byte repeats) IncAmplitude Gets the inc amplitude command.	12
protected abstract getIncFrequencyCommand (byte groupNum) IncFrequency Gets the inc frequency command.	12
protected abstract getIncPulseWidthCommand (byte progNum) IncPulseWidth Gets the inc pulse width command.	12
protected abstract getResetCycleCommand () ResetCycle Gets the reset cycle command.	12
protected abstract getRestoreClinicianSettingsCommand (byte groupNumber) RestoreClinicianSettings Gets the restore clinician settings command.	13
protected abstract getSendTriggerCommand () SendTrigger Gets the send trigger command.	13
protected abstract getSetActiveGroupCommand (byte gpNum) SetActiveGroup Gets the sets the active group command.	13
protected abstract getSetNexusConfigurationCommand (byte maintSessTimeoutSec, byte hostSessTimeoutMin) SetNexusConfiguration Gets the sets the nexus configuration command.	13
protected abstract getStartRealTimeCommand (boolean channel1) StartRealTime Gets the start real time command.	13

protected abstract <i>StartSensing</i>	getStartSensingCommand() Gets the start sensing command.	14
protected abstract <i>StopRealTime</i>	getStopRealTimeCommand() Gets the stop real time command.	14
protected abstract <i>StopSensing</i>	getStopSensingCommand() Gets the stop sensing command.	14
protected abstract <i>TherapyOff</i>	getTherapyOffCommand() Gets the therapy off command.	14
protected abstract <i>TherapyOn</i>	getTherapyOnCommand() Gets the therapy on command.	14
static void	setProvider (<i>CommandAccessor</i> provider) Sets the provider.	10

Constructor Detail

CommandAccessor

```
public CommandAccessor()
```

Method Detail

setProvider

```
public static void setProvider(CommandAccessor provider)
```

Sets the provider.

Parameters:

provider - the new provider

getCycleStimCommand

```
protected abstract CycleStim getCycleStimCommand()
```

Gets the cycle stim command.

Returns:

the cycle stim command

getDecAmplitudeCommand

```
protected abstract DecAmplitude getDecAmplitudeCommand(byte progNum,  
                                                         byte numSteps,  
                                                         byte repeats)
```

Gets the dec amplitude command.

Parameters:

progNum - the prog num

numSteps - the num steps

repeats - the repeats
Returns:
the dec amplitude command

getDecFrequencyCommand

protected abstract [DecFrequency](#) getDecFrequencyCommand(byte groupNum)

Gets the dec frequency command.

Parameters:

groupNum - the group num

Returns:
the dec frequency command

getDecPulseWidthCommand

protected abstract [DecPulseWidth](#) getDecPulseWidthCommand(byte progNum)

Gets the dec pulse width command.

Parameters:

progNum - the prog num

Returns:
the dec pulse width command

getGetInsInfoCommand

protected abstract [GetInsInfo](#) getGetInsInfoCommand()

Gets the gets the ins info command.

Returns:
the gets the ins info command

getGetNexusStatusCommand

protected abstract [GetNexusStatus](#) getGetNexusStatusCommand()

Gets the gets the nexus status command.

Returns:
the gets the nexus status command

getGetRealTimeDataCommand

protected abstract [GetRealtimeData](#) getGetRealTimeDataCommand()

Gets the gets the real time data command.

Returns:
the gets the real time data command

getIncAmplitudeCommand

```
protected abstract IncAmplitude getIncAmplitudeCommand(byte progNum,  
                                                         byte numSteps,  
                                                         byte repeats)
```

Gets the inc amplitude command.

Parameters:
progNum - the prog num
numSteps - the num steps
repeats - the repeats

Returns:
the inc amplitude command

getIncFrequencyCommand

```
protected abstract IncFrequency getIncFrequencyCommand(byte groupNum)
```

Gets the inc frequency command.

Parameters:
groupNum - the group num

Returns:
the inc frequency command

getIncPulseWidthCommand

```
protected abstract IncPulseWidth getIncPulseWidthCommand(byte progNum)
```

Gets the inc pulse width command.

Parameters:
progNum - the prog num

Returns:
the inc pulse width command

getResetCycleCommand

```
protected abstract ResetCycle getResetCycleCommand()
```

Gets the reset cycle command.

Returns:
the reset cycle command

getRestoreClinicianSettingsCommand

```
protected abstract RestoreClinicianSettings getRestoreClinicianSettingsCommand(byte groupNumber)
```

Gets the restore clinician settings command.

Parameters:

groupNumber - the group number

Returns:

the restore clinician settings command

getSendTriggerCommand

```
protected abstract SendTrigger getSendTriggerCommand()
```

Gets the send trigger command.

Returns:

the send trigger command

getSetActiveGroupCommand

```
protected abstract SetActiveGroup getSetActiveGroupCommand(byte gpNum)
```

Gets the sets the active group command.

Parameters:

gpNum - the gp num

Returns:

the sets the active group command

getSetNexusConfigurationCommand

```
protected abstract SetNexusConfiguration getSetNexusConfigurationCommand(byte maintSessTimeoutSec,
                                                                              byte hostSessTimeoutMin)
in)
```

Gets the sets the nexus configuration command.

Parameters:

maintSessTimeoutSec - the maint sess timeout sec

hostSessTimeoutMin - the host sess timeout min

Returns:

the sets the nexus configuration command

getStartRealTimeCommand

```
protected abstract StartRealTime getStartRealTimeCommand(boolean channel1)
```

Gets the start real time command.

Parameters:

channel1 - the channel1

Returns:

the start real time command

getStartSensingCommand

```
protected abstract StartSensing getStartSensingCommand()
```

Gets the start sensing command.

Returns:

the start sensing command

getStopRealTimeCommand

```
protected abstract StopRealTime getStopRealTimeCommand()
```

Gets the stop real time command.

Returns:

the stop real time command

getStopSensingCommand

```
protected abstract StopSensing getStopSensingCommand()
```

Gets the stop sensing command.

Returns:

the stop sensing command

getTherapyOffCommand

```
protected abstract TherapyOff getTherapyOffCommand()
```

Gets the therapy off command.

Returns:

the therapy off command

getTherapyOnCommand

```
protected abstract TherapyOn getTherapyOnCommand()
```

Gets the therapy on command.

Returns:

the therapy on command

Enum ConnectionStatus

[mdt.neuro.nexus](#)

```
java.lang.Object
├─ java.lang.Enum<ConnectionStatus>
│   └─ mdt.neuro.nexus.ConnectionStatus
```

All Implemented Interfaces:

[Comparable](#)<[ConnectionStatus](#)>, [Serializable](#)

```
public enum ConnectionStatus
extends Enum<ConnectionStatus>
```

The Enum ConnectionStatus.

Enum Constant Summary	Page
Connected Communication established.	15
Connecting In the process of establishing communication.	15
Disconnected Communication is not established.	15

Method Summary	Page
<small>static ConnectionStatus</small> valueOf (String name)	16
<small>static ConnectionStatus[]</small> values ()	16

Enum Constant Detail

Connected

```
public static final ConnectionStatus Connected
```

Communication established.

Connecting

```
public static final ConnectionStatus Connecting
```

In the process of establishing communication.

Disconnected

```
public static final ConnectionStatus Disconnected
```

Communication is not established.

Method Detail

values

```
public static ConnectionStatus[] values()
```

valueOf

```
public static ConnectionStatus valueOf(String name)
```

Class `DataPacketAccessor`

`mdt.neuro.nexus``java.lang.Object``└ mdt.neuro.nexus.DataPacketAccessor`

```
abstract public class DataPacketAccessor
extends Object
```

The Class `DataPacketAccessor`. This class defines methods to allow access to `DataPacket` methods to other internal classes of the API

Constructor Summary	Page
<code>DataPacketAccessor</code> ()	17

Method Summary		Page
protected abstract void	setNumMissedPatterns (DataPacket pkt, short misses) Sets the num missed packets.	18
protected abstract void	setPacketNum (DataPacket pkt, int index, short num) Sets the packet num.	17
static void	setProvider (DataPacketAccessor provider) Sets the provider.	17

Constructor Detail

`DataPacketAccessor`

```
public DataPacketAccessor()
```

Method Detail

`setProvider`

```
public static void setProvider(DataPacketAccessor provider)
```

Sets the provider.

Parameters:

provider - the new provider

`setPacketNum`

```
protected abstract void setPacketNum(DataPacket pkt,
                                       int index,
                                       short num)
```

Sets the packet num.

Parameters:

pkt - the pkt
index - the index
num - the num

setNumMissedPatterns

```
protected abstract void setNumMissedPatterns(DataPacket pkt,  
                                              short misses)
```

Sets the num missed packets.

Parameters:

pkt - the pkt
misses - the misses

Interface IDisposable

[mdt.neuro.nexus](#)

All Known Implementing Classes:

[ObservableConnection](#), [SerialConnection](#)

```
public interface IDisposable
```

The Interface IDisposable.

Method Summary		Page
void	dispose() Dispose of this instance.	19

Method Detail

dispose

```
void dispose()
```

Dispose of this instance.

Class NexusInstrument

mdt.neuro.nexus

```

java.lang.Object
├── java.util.Observable
│   └── mdt.neuro.nexus.NexusInstrument

```

All Implemented Interfaces:

Observer

```

public class NexusInstrument
extends Observable
implements Observer

```

The Class NexusInstrument. This class is the main abstraction of the Nexus-2 system interface to the host application.

Method Summary		Page
List<Byte>	arbitraryCommand (List<Byte> bytes, int timeoutMs) Sends a command of arbitrary bytes.	23
double[][]	computeTemplate (DataPacket[] pkts) Computes a template based on passed in data packets that can then be subtracted from other data packets to reduce artifact induced by telemetry.	29
int	connect (ObservableConnection conn) Connect this instance and initializes the internal API classes required for communication.	22
int	cycleStim () Cycle stim on the device (off then back on).	28
double	decStimParameter (TherapyParameter type, byte progNum) Decrements the given stimulation parameter.	24
double	decStimParameter (TherapyParameter type, byte progNum, byte numSteps, byte numRepeats) Decrements the given stimulation parameter.	24
boolean	disconnect () Disconnects this instance.	22
void	dispose () Disposes of this instance.	22
DataPacket	getDataPacket () Gets a data packet from the INS.	25
InsInfo	getInsInfo () Gets the INS info, including stimulation settings and sensing status.	25
static NexusInstrument	getInstance () Gets the single instance of NexusInstrument.	22
int	getLastInsResponseCode () Gets the last ins response code.	22
ResponseCodes	getLastNexusResponseCode () Gets the last nexus response code.	22
NexusStatus	getNexusStatus () Gets the status of the Nexus-D system.	25

String	getVersion() Gets the version of the API.	21
double	incStimParameter (TherapyParameter type, byte progNum) Increments the given stimulation parameter.	23
double	incStimParameter (TherapyParameter type, byte progNum, byte numSteps, byte numRepeats) Increments the given stimulation parameter.	24
double	pulseStimParameter (byte progNum, byte numSteps, byte numDelay) Sends a batched inc/dec to allow for fine control of the number of stimulation pulses that are delivered at a given level.	29
int	resetCycle() Resets the stimulation cycle.	28
int	restoreClinicianSettings (byte groupNumber) Restores the clinician settings saved in the INS.	28
int	sendTrigger() Send a trigger to the INS.	28
int	setActiveGroup (byte groupNum) Sets the active group.	27
int	setNexusConfiguration (byte maintSessionTimeoutSec, byte hostSessionTimeoutMin) Sets the nexus configuration.	25
int	startDataSession() Starts a data session with the INS.	26
int	startDataSession (boolean channelOne) Starts a data session with the INS.	26
int	startSensing() Enables the sensing feature in the INS.	26
int	stopDataSession() Stops a data session.	27
int	stopSensing() Disables the sensing feature in the INS.	26
int	therapyOff() Turns therapy off.	27
int	therapyOn() Turns therapy on.	27
String	toString()	23
void	update (Observable o, Object arg)	23

Method Detail

getVersion

```
public String getVersion()
```

Gets the version of the API.

Returns:

the version

getLastNexusResponseCode

```
public ResponseCodes getLastNexusResponseCode()
```

Gets the last nexus response code.

Returns:
the nexusResponseCode

getLastInsResponseCode

```
public int getLastInsResponseCode()
```

Gets the last ins response code.

Returns:
the insResponseCode

getInstance

```
public static synchronized NexusInstrument getInstance()
```

Gets the single instance of NexusInstrument.

Returns:
single instance of NexusInstrument

connect

```
public int connect(ObservableConnection conn)
```

Connect this instance and initializes the internal API classes required for communication.

Parameters:
`conn` - the observable connection to connect to (e.g. serial connection)

Returns:
result of the operation (0 for success)

disconnect

```
public boolean disconnect()
```

Disconnects this instance.

Returns:
true, if disconnect was successful

dispose

```
public void dispose()
```

Disposes of this instance.

update

```
public void update(Observable o,  
                  Object arg)
```

Specified by:

update in interface *Observer*

toString

```
public String toString()
```

Overrides:

toString in class *Object*

arbitraryCommand

```
public List<Byte> arbitraryCommand(List<Byte> bytes,  
                                   int timeoutMs)  
    throws InterruptedException
```

Sends a command of arbitrary bytes.

Parameters:

bytes - the bytes that compose the command. This includes the command code and parameter bytes but does not include the CRC.
timeoutMs - the timeout ms

Returns:

the object returned

Throws:

InterruptedException - the interrupted exception

incStimParameter

```
public double incStimParameter(TherapyParameter type,  
                               byte progNum)  
    throws InterruptedException
```

Increments the given stimulation parameter.

Parameters:

type - the parameter type
progNum - the program number

Returns:

the resulting *TherapyParameter* value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known *AmplitudeResolution* (V or CC). If this is changed without sending a *GetInsInfo* command to give the API knowledge of the change, the returned value will be incorrect

Throws:

InterruptedException - the interrupted exception

incStimParameter

```
public double incStimParameter(TherapyParameter type,
                               byte progNum,
                               byte numSteps,
                               byte numRepeats)
    throws InterruptedException
```

Increments the given stimulation parameter.

Parameters:

type - the parameter type
progNum - the program number
numSteps - the number of steps to increment (only for amplitude)
numRepeats - the num repeats (only for amplitude)

Returns:

the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Throws:

[InterruptedException](#) - the interrupted exception

decStimParameter

```
public double decStimParameter(TherapyParameter type,
                               byte progNum)
    throws InterruptedException
```

Decrements the given stimulation parameter.

Parameters:

type - the parameter type
progNum - the program number

Returns:

the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Throws:

[InterruptedException](#) - the interrupted exception

decStimParameter

```
public double decStimParameter(TherapyParameter type,
                               byte progNum,
                               byte numSteps,
                               byte numRepeats)
    throws InterruptedException
```

Decrements the given stimulation parameter.

Parameters:

type - the parameter type
progNum - the program number
numSteps - the number of steps to decrement (only for amplitude)

numRepeats - the num repeats (only for amplitude)

Returns:

the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect.

Throws:

[InterruptedException](#) - the interrupted exception

getNexusStatus

```
public NexusStatus getNexusStatus ()  
    throws InterruptedException
```

Gets the status of the Nexus-D system.

Returns:

the Nexus status object or null in the case of an error

Throws:

[InterruptedException](#) - the interrupted exception

getInsInfo

```
public InsInfo getInsInfo ()  
    throws InterruptedException
```

Gets the INS info, including stimulation settings and sensing status.

Returns:

the INS info object or null in the case of an error

Throws:

[InterruptedException](#) - the interrupted exception

getDataPacket

```
public DataPacket getDataPacket ()  
    throws InterruptedException
```

Gets a data packet from the INS. The contents of the data packet are dependent upon the sensing settings configured in the INS. These sensing settings cannot be configured by the Nexus-D system.

Returns:

the data packet or null in the case of an error

Throws:

[InterruptedException](#) - the interrupted exception

setNexusConfiguration

```
public int setNexusConfiguration (byte maintSessionTimeoutSec,  
    byte hostSessionTimeoutMin)  
    throws InterruptedException
```

Sets the nexus configuration.

Parameters:

maintSessionTimeoutSec - the maintenance session timeout in seconds
hostSessionTimeoutMin - the host session timeout minutes

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

startSensing

```
public int startSensing()  
    throws InterruptedException
```

Enables the sensing feature in the INS.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

stopSensing

```
public int stopSensing()  
    throws InterruptedException
```

Disables the sensing feature in the INS.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

startDataSession

```
public int startDataSession(boolean channelOne)  
    throws InterruptedException
```

Starts a data session with the INS.

Parameters:

channelOne - If 422 Hz is the TD sampling frequency, channel 1 will be the channel used based on this boolean

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

startDataSession

```
public int startDataSession()  
    throws InterruptedException
```

Starts a data session with the INS. If 422 Hz is the TD sampling frequency, channel 1 will be the default channel used

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

`InterruptedException` - the interrupted exception

stopDataSession

```
public int stopDataSession()  
    throws InterruptedException
```

Stops a data session.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

`InterruptedException` - the interrupted exception

therapyOn

```
public int therapyOn()  
    throws InterruptedException
```

Turns therapy on.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

`InterruptedException` - the interrupted exception

therapyOff

```
public int therapyOff()  
    throws InterruptedException
```

Turns therapy off.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

`InterruptedException` - the interrupted exception

setActiveGroup

```
public int setActiveGroup(byte groupNum)  
    throws InterruptedException
```

Sets the active group.

Parameters:

groupNum - the group number to activate

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

restoreClinicianSettings

```
public int restoreClinicianSettings(byte groupNumber)
    throws InterruptedException
```

Restores the clinician settings saved in the INS.

Parameters:

groupNumber - the group number to restore (00 for all)

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

sendTrigger

```
public int sendTrigger()
    throws InterruptedException
```

Send a trigger to the INS.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

cycleStim

```
public int cycleStim()
    throws InterruptedException
```

Cycle stim on the device (off then back on).

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

resetCycle

```
public int resetCycle()
    throws InterruptedException
```

Resets the stimulation cycle.

Returns:

the result of the operation (INS response code) or -1 for a timeout

Throws:

[InterruptedException](#) - the interrupted exception

pulseStimParameter

```
public double pulseStimParameter(byte progNum,  
                                byte numSteps,  
                                byte numDelay)  
    throws InterruptedException
```

Sends a batched inc/dec to allow for fine control of the number of stimulation pulses that are delivered at a given level.

Parameters:

progNum - the prog num
numSteps - the num steps
numDelay - the num delay

Returns:

the new amplitude

Throws:

[InterruptedException](#) - the interrupted exception

computeTemplate

```
public double[][] computeTemplate(DataPacket[] pkts)
```

Computes a template based on passed in data packets that can then be subtracted from other data packets to reduce artifact induced by telemetry. At least 20 packets must be passed in to compute a valid template.

Parameters:

pkts - the pkts

Returns:

the double[][] template

Class ObservableConnection

mdt.neuro.nexus

```

java.lang.Object
├── java.util.Observable
│   └── mdt.neuro.nexus.ObservableConnection

```

All Implemented Interfaces:

IDisposable

Direct Known Subclasses:

SerialConnection

```

abstract public class ObservableConnection
extends Observable
implements IDisposable

```

The Class ObservableConnection.

Field Summary		Page
protected ConnectionStatus	status The status.	31

Constructor Summary		Page
ObservableConnection ()		31

Method Summary		Page
protected abstract int	connect () Connect.	31
protected abstract boolean	disconnect () Disconnect.	31
void	dispose () Disposes of this instance.	32
protected abstract int	getInputBufferBytesCount () Gets the input buffer bytes count.	32
ConnectionStatus	getStatus () Gets the status of this connection.	31
protected abstract byte[]	read () Read from the connection.	32
protected void	receive () Receives data from the connection.	32
protected abstract int	write (byte[] bufferedData, int count) Writes data to the connection.	32
protected int	writeBytes (byte[] bufferedData, int count) Write to the connection.	31

Field Detail

status

protected `ConnectionStatus` **status**

The status.

Constructor Detail

ObservableConnection

public `ObservableConnection()`

Method Detail

getStatus

public `ConnectionStatus` **getStatus()**

Gets the status of this connection.

Returns:
the status

connect

protected abstract int **connect()**

Connect.

Returns:
0, if successful. Error code otherwise

disconnect

protected abstract boolean **disconnect()**

Disconnect.

Returns:
true, if successful

writeBytes

protected int **writeBytes**(byte[] bufferedData,
int count)

Write to the connection.

Parameters:

bufferedData - the buffered data
count - the count

Returns:

the int

write

```
protected abstract int write(byte[] bufferedData,  
                             int count)
```

Writes data to the connection.

Parameters:

bufferedData - the buffered data
count - the count

Returns:

the int

receive

```
protected void receive()
```

Receives data from the connection. Processes the data according to the Nexus-D protocol for forming header and payload

read

```
protected abstract byte[] read()
```

Read from the connection.

Returns:

the byte[] bytes that were read

getInputBufferBytesCount

```
protected abstract int getInputBufferBytesCount()
```

Gets the input buffer bytes count.

Returns:

the input buffer bytes count

dispose

```
public void dispose()
```

Disposes of this instance.

Specified by:

`dispose` in interface `IDisposable`

Class ParameterizedRunnable

mdt.neuro.nexus

java.lang.Object
└ mdt.neuro.nexus.ParameterizedRunnable

All Implemented Interfaces:
Runnable

abstract public class **ParameterizedRunnable**
extends [Object](#)
implements [Runnable](#)

The Class ParameterizedRunnable.

Constructor Summary		Page
ParameterizedRunnable (Object parameter)	Instantiates a new parameterized runnable.	34

Method Summary		Page
void	run ()	34
abstract void	runWithParameter (Object param) Executes with the given parameter.	34

Constructor Detail

ParameterizedRunnable

public **ParameterizedRunnable**([Object](#) parameter)

Instantiates a new parameterized runnable.

Parameters:
parameter - the parameter

Method Detail

runWithParameter

public abstract void **runWithParameter**([Object](#) param)

Executes with the given parameter.

Parameters:
param - the parameter

run

public void **run**()

Specified by:

`run` in interface `Runnable`

Enum ResponseCodes

[mdt.neuro.nexus](#)

```
java.lang.Object
└─ java.lang.Enum<ResponseCodes>
    └─ mdt.neuro.nexus.ResponseCodes
```

All Implemented Interfaces:

[Comparable<ResponseCodes>](#), [Serializable](#)

```
public enum ResponseCodes
extends Enum<ResponseCodes>
```

The Enum ResponseCodes.

Enum Constant Summary	Page
BATT_DEPLETED STS2 battery has been depleted.	38
HEADER_CRC_ERROR Header crc error.	37
INS_POR_INDICATED INS POR indicated.	38
INVALID_DATA	38
INVALID_FRAME_ID Invalid frame id.	37
INVALID_FRAME_TYPE Invalid frame type.	37
INVALID_PAYLOAD_LENGTH Invalid payload length.	37
MESSAGE_INCOMPLETE Message Incomplete.	37
NO_RESPONSE No response received from the STS2.	38
NOT_CONNECTED No connection to the STS2.	38
PAYLOAD_CRC_ERROR Payload crc error.	37
PREV_CMD_BUSY Previous command busy.	38
SUCCESS Success.	37

Method Summary	Page
static ResponseCodes valueOf (String name)	38
static ResponseCodes[] values ()	38

Enum Constant Detail

SUCCESS

```
public static final ResponseCodes SUCCESS
```

Success.

PAYLOAD_CRC_ERROR

```
public static final ResponseCodes PAYLOAD_CRC_ERROR
```

Payload crc error.

INVALID_FRAME_TYPE

```
public static final ResponseCodes INVALID_FRAME_TYPE
```

Invalid frame type.

MESSAGE_INCOMPLETE

```
public static final ResponseCodes MESSAGE_INCOMPLETE
```

Message Incomplete.

INVALID_FRAME_ID

```
public static final ResponseCodes INVALID_FRAME_ID
```

Invalid frame id.

INVALID_PAYLOAD_LENGTH

```
public static final ResponseCodes INVALID_PAYLOAD_LENGTH
```

Invalid payload length.

HEADER_CRC_ERROR

```
public static final ResponseCodes HEADER_CRC_ERROR
```

Header crc error.

PREV_CMD_BUSY

```
public static final ResponseCodes PREV_CMD_BUSY
```

Previous command busy.

INS_POR_INDICATED

```
public static final ResponseCodes INS_POR_INDICATED
```

INS POR indicated.

BATT_DEPLETED

```
public static final ResponseCodes BATT_DEPLETED
```

STS2 battery has been depleted.

INVALID_DATA

```
public static final ResponseCodes INVALID_DATA
```

NO_RESPONSE

```
public static final ResponseCodes NO_RESPONSE
```

No response received from the STS2.

NOT_CONNECTED

```
public static final ResponseCodes NOT_CONNECTED
```

No connection to the STS2.

Method Detail

values

```
public static ResponseCodes[] values()
```

valueOf

```
public static ResponseCodes valueOf(String name)
```

Class SerialConnection

mdt.neuro.nexus

```

java.lang.Object
├── java.util.Observable
│   └── mdt.neuro.nexus.ObservableConnection
│       └── mdt.neuro.nexus.SerialConnection

```

All Implemented Interfaces:

IDisposable

```

public class SerialConnection
extends ObservableConnection

```

The Class Connection.

Fields inherited from class mdt.neuro.nexus.ObservableConnection

status

Constructor Summary

SerialConnection(String portName)
Instantiates a new connection.

Page

39

Method Summary

		Page
protected int	connect () Connects this serial instance to the given port.	40
protected boolean	disconnect () Disconnects this instance.	40
protected int	getInputBufferBytesCount () Gets the input buffer bytes count.	41
protected byte[]	read () Read from the connection.	41
void	serialEvent (SerialPortEvent event)	41
boolean	setPortName (String portName) Sets the port name and initializes the port if it does not exist.	40
protected int	write (byte[] toSend, int count) Writes data to the connection.	40

Methods inherited from class mdt.neuro.nexus.ObservableConnection

dispose, getStatus, receive, writeBytes

Constructor Detail

SerialConnection

```
public SerialConnection(String portName)
```

Instantiates a new connection.

Parameters:

portName - the port name

Method Detail

setPortName

```
public boolean setPortName(String portName)
    throws SerialPortException
```

Sets the port name and initializes the port if it does not exist. If a port is already open, it will be closed and a new one will be initialized using the new name

Parameters:

portName - the port name

Returns:

true, if successful

Throws:

SerialPortException - the serial port exception

connect

```
protected int connect()
```

Connects this serial instance to the given port.

Overrides:

[connect](#) in class [ObservableConnection](#)

Returns:

true, if successful

disconnect

```
protected boolean disconnect()
```

Disconnects this instance.

Overrides:

[disconnect](#) in class [ObservableConnection](#)

Returns:

true, if successful

write

```
protected int write(byte[] toSend,
    int count)
```

Writes data to the connection.

Overrides:

[write](#) in class [ObservableConnection](#)

Parameters:

count - the count

Returns:
the int

serialEvent

public void **serialEvent**(SerialPortEvent event)

read

protected byte[] **read**()

Read from the connection.

Overrides:
[read](#) in class [ObservableConnection](#)

Returns:
the byte[] bytes that were read

getInputBufferBytesCount

protected int **getInputBufferBytesCount**()

Gets the input buffer bytes count.

Overrides:
[getInputBufferBytesCount](#) in class [ObservableConnection](#)

Returns:
the input buffer bytes count

Class SerialPortErrors

mdt.neuro.nexus

java.lang.Object
└ mdt.neuro.nexus.SerialPortErrors

public class **SerialPortErrors**
extends **Object**

The Class SerialPortErrors.

Field Summary		Page
static int	PORT_BUSY The port was busy.	42
static int	PORT_NOT_FOUND The port was not found.	42
static int	PORT_NULL The port was not defined.	42
static int	UNKNOWN The unknown.	42

Field Detail

PORT_NOT_FOUND

public static final int **PORT_NOT_FOUND**

The port was not found.

PORT_BUSY

public static final int **PORT_BUSY**

The port was busy.

PORT_NULL

public static final int **PORT_NULL**

The port was not defined.

UNKNOWN

public static final int **UNKNOWN**

The unknown.

Class ThreadedNexusInstrument

mdt.neuro.nexus

java.lang.Object

└ mdt.neuro.nexus.ThreadedNexusInstrument

```
public class ThreadedNexusInstrument
extends Object
```

The Class ThreadedNexusInstrument.

Constructor Summary	Page
ThreadedNexusInstrument()	44

Method Summary	Page
void arbitraryCommand (List<Byte> bytesToSend, int timeoutMs) Sends an command of arbitrary bytes.	48
void cycleStim () Cycles the stimulation.	48
void decStimParameter (TherapyParameter type, byte progNum) Decrement the given stimulation parameter.	45
void decStimParameter (TherapyParameter type, byte progNum, byte numSteps, byte numRepeats) Decrement the given stimulation parameter.	46
void getDataPacket () Gets a data packet from the INS.	46
void getInsInfo () Gets the INS info, including stimulation settings and sensing status.	46
int getLastInsResponseCode () Gets the last ins response code.	45
ResponseCodes getLastNexusResponseCode () Gets the last nexus response code.	44
void getNexusStatus () Gets the status of the Nexus-D system.	46
Object getThreadSafeReturnVal () Gets the thread safe return val.	44
void incStimParameter (TherapyParameter type, byte progNum) Increment the given stimulation parameter.	45
void incStimParameter (TherapyParameter type, byte progNum, byte numSteps, byte numRepeats) Increment the given stimulation parameter.	45
void pulseStimParameter (byte progNum, byte numSteps, byte numDelay) Sends a batched inc/dec to allow for fine control of the number of stimulation pulses that are delivered at a given level.	49
void resetCycle () Reset the stimulation cycle.	49
void restoreClinicianSettings (byte group) Restores the clinician settings saved in the INS.	48

void	sendTrigger () Send a trigger to the INS.	48
void	setActiveGroup (byte group) Sets the active group.	48
void	setNexusConfiguration (byte maintSessionTimeoutSec, byte hostSessionTimeoutMin) Sets the nexus configuration.	46
void	startDataSession () Starts a data session with the INS.	47
void	startDataSession (boolean channelOne) Starts a data session with the INS.	47
void	startSensing () Enables the sensing feature in the INS.	47
void	stopDataSession () Stops a data session.	47
void	stopSensing () Disables the sensing feature in the INS.	47
void	therapyOff () Turns therapy off.	48
void	therapyOn () Turns therapy on.	47

Constructor Detail

ThreadedNexusInstrument

```
public ThreadedNexusInstrument ()
```

Method Detail

getThreadSafeReturnVal

```
public Object getThreadSafeReturnVal ()
```

Gets the thread safe return val.

Returns:
the thread safe return val

getLastNexusResponseCode

```
public ResponseCodes getLastNexusResponseCode ()
```

Gets the last nexus response code.

Returns:
the nexusResponseCode

getLastInsResponseCode

```
public int getLastInsResponseCode()
```

Gets the last ins response code.

Returns:
the insResponseCode

incStimParameter

```
public void incStimParameter(TherapyParameter type,  
                             byte progNum)
```

Increment the given stimulation parameter. Sets the returned value to the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Parameters:
type - the type
progNum - the prog num

incStimParameter

```
public void incStimParameter(TherapyParameter type,  
                             byte progNum,  
                             byte numSteps,  
                             byte numRepeats)
```

Increment the given stimulation parameter. Sets the returned value to the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Parameters:
type - the type
progNum - the prog num
numSteps - the num steps
numRepeats - the num repeats

decStimParameter

```
public void decStimParameter(TherapyParameter type,  
                             byte progNum)
```

Decrement the given stimulation parameter. Sets the returned value to the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Parameters:
type - the type
progNum - the prog num

decStimParameter

```
public void decStimParameter(TherapyParameter type,
                             byte progNum,
                             byte numSteps,
                             byte numRepeats)
```

Decrement the given stimulation parameter. Sets the returned value to the resulting [TherapyParameter](#) value as a result of the operation. -1 if the operation failed. The resulting amplitude value is calculated based on the last known [AmplitudeResolution](#) (V or CC). If this is changed without sending a [GetInsInfo](#) command to give the API knowledge of the change, the returned value will be incorrect

Parameters:

type - the type
progNum - the prog num
numSteps - the num steps
numRepeats - the num repeats

getNexusStatus

```
public void getNexusStatus()
```

Gets the status of the Nexus-D system. Sets the returned value to the nexus status object or null in the case of an error

getInsInfo

```
public void getInsInfo()
```

Gets the INS info, including stimulation settings and sensing status. Sets the returned value to the ins info object or null in the case of an error

getDataPacket

```
public void getDataPacket()
```

Gets a data packet from the INS. The contents of the data packet are dependent upon the sensing settings configured in the INS. These sensing settings cannot be configured by the Nexus-D system. Sets the returned value to the data packet or null in the case of an error

setNexusConfiguration

```
public void setNexusConfiguration(byte maintSessionTimeoutSec,
                                   byte hostSessionTimeoutMin)
```

Sets the nexus configuration. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

Parameters:

maintSessionTimeoutSec - the maint session timeout sec
hostSessionTimeoutMin - the host session timeout min

startSensing

```
public void startSensing()
```

Enables the sensing feature in the INS. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

stopSensing

```
public void stopSensing()
```

Disables the sensing feature in the INS. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

startDataSession

```
public void startDataSession()
```

Starts a data session with the INS. If 422 Hz is the TD sampling frequency, channel 1 will be the default channel used. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

startDataSession

```
public void startDataSession(boolean channelOne)
```

Starts a data session with the INS. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

Parameters:

`channelOne` - If 422 Hz is the TD sampling frequency, channel 1 will be the channel used based on this boolean

stopDataSession

```
public void stopDataSession()
```

Stops a data session. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

therapyOn

```
public void therapyOn()
```

Turns therapy on. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

therapyOff

```
public void therapyOff()
```

Turns therapy off. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

setActiveGroup

```
public void setActiveGroup(byte group)
```

Sets the active group. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

Parameters:

group - the new active group

restoreClinicianSettings

```
public void restoreClinicianSettings(byte group)
```

Restores the clinician settings saved in the INS. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

Parameters:

group - the group

sendTrigger

```
public void sendTrigger()
```

Send a trigger to the INS. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

cycleStim

```
public void cycleStim()
```

Cycles the stimulation. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

arbitraryCommand

```
public void arbitraryCommand(List<Byte> bytesToSend,  
                             int timeoutMs)
```

Sends an command of arbitrary bytes. Sets the returned value to the returned bytes from the operation

Parameters:

bytesToSend - the bytes to send

timeoutMs - the timeout ms

resetCycle

```
public void resetCycle()
```

Reset the stimulation cycle. Sets the returned value to the result of the operation (INS response code) or -1 for a timeout

pulseStimParameter

```
public void pulseStimParameter(byte progNum,  
                                byte numSteps,  
                                byte numDelay)
```

Sends a batched inc/dec to allow for fine control of the number of stimulation pulses that are delivered at a given level.

Parameters:

progNum - the prog num
numSteps - the num steps
numDelay - the num delay

Class ThreadManager

mdt.neuro.nexus

java.lang.Object

└ mdt.neuro.nexus.ThreadManager

All Implemented Interfaces:

Observer

```
public class ThreadManager
extends Object
implements Observer
```

The Class ThreadManager. This class is a singleton class that allows a Runnable to be executed on a worker thread. It supports the following functionality: Monitor its execution and abort an operation in progress. This method interrupts the Runnable that is being executed. The Runnable should handle the InterruptedException that will be thrown gracefully to ensure consistent operation. The ThreadManager only allows one operation to be executed at a time, and does not support the queuing of operations.

Method Summary		Page
void	abortOperation() Abort operation.	51
boolean	doOperation (Runnable toExecute) Does operation specified by the Runnable on the worker thread if it is available.	50
static ThreadManager	getInstance() Gets the single instance of ThreadManager.	50
boolean	isExecuting()	50
void	update (Observable o, Object arg)	51

Method Detail

isExecuting

```
public boolean isExecuting()
```

getInstance

```
public static ThreadManager getInstance()
```

Gets the single instance of ThreadManager.

Returns:

single instance of ThreadManager

doOperation

```
public boolean doOperation(Runnable toExecute)
```

Does operation specified by the Runnable on the worker thread if it is available.

Parameters:

toExecute - the Runnable to execute

Returns:

true, if the operation was started on the worker thread

abortOperation

```
public void abortOperation()
```

Abort operation. This will interrupt the worker thread. The implemented runnable must handle a potential `InterruptedException` gracefully

update

```
public synchronized void update(Observable o,  
                                Object arg)
```

Specified by:

update in interface `Observer`

Class Util

[mdt.neuro.nexus](#)

```
java.lang.Object
└─ mdt.neuro.nexus.Util
```

```
public class Util
extends Object
```

The Class Util.

Constructor Summary	Page
Util()	52

Method Summary		Page
static String	toHex (byte[] array) Converts a byte array to a hex string.	52
static String	toHex (Byte[] array) Converts a byte array to a hex string.	53
static Byte[]	toObject (byte[] array) Converts a byte array to a Byte array.	52

Constructor Detail

Util

```
public Util()
```

Method Detail

toObject

```
public static Byte[] toObject(byte[] array)
```

Converts a byte array to a Byte array.

Parameters:

array - the array

Returns:

the byte[]

toHex

```
public static String toHex(byte[] array)
```

Converts a byte array to a hex string.

Parameters:

array - the array

Returns:
the string

toHex

```
public static String toHex(Byte[] array)
```

Converts a byte array to a hex string.

Parameters:
array - the array

Returns:
the string

Package mdt.neuro.nexus.commands

Class Summary		Page
CommandCodes	The Class CommandCodes.	54
CommandResponseCodes	The Class CommandResponseCodes.	60
CycleStim	The Class CycleStim.	69
DataStructureAccessor	The Class DataStructureAccessor.	70
DecAmplitude	The Class DecAmplitude.	72
DecFrequency	The Class DecFrequency.	74
DecPulseWidth	The Class DecPulseWidth.	76
GetInsInfo	The Class GetInsInfo.	78
GetNexusStatus	The Class GetNexusStatus.	80
GetRealtimeData	The Class GetRealtimeData.	82
IncAmplitude	The Class IncAmplitude.	84
IncFrequency	The Class IncFrequency.	86
IncPulseWidth	The Class IncPulseWidth.	88
PulseStim	The Class PulseStim.	90
ResetCycle		92
RestoreClinicianSettings	The Class RestoreClinicianSettings.	93
SendTrigger	The Class SendTrigger.	94
SetActiveGroup	The Class SetActiveGroup.	95
SetNexusConfiguration	The Class SetNexusConfiguration.	96
StartRealTime	The Class StartRealTime.	97
StartSensing	The Class StartSensing.	98
StopRealTime	The Class StopRealTime.	99
StopSensing	The Class StopSensing.	100
TherapyOff	The Class TherapyOff.	101
TherapyOn	The Class TherapyOn.	102

Class CommandCodes

[mdt.neuro.nexus.commands](#)

[java.lang.Object](#)

└ [mdt.neuro.nexus.commands.CommandCodes](#)

```
public class CommandCodes
extends Object
```

The Class CommandCodes.

Field Summary		Page
static short	CYCLE_STIM The cycle stim.	57
static short	DEC_AMP The dec amp.	58
static short	DEC_FREQ The dec freq.	58
static short	DEC_PULSEWIDTH The dec pulsewidth.	58
static short	GET_INS_INFO The get ins info.	56
static short	GET_REALTIME_DATA The get realtime data.	57
static short	GET_STATUS The get status.	57
static short	INC_AMP The inc amp.	58
static short	INC_FREQ The inc freq.	58
static short	INC_PULSEWIDTH The inc pulsewidth.	58
static short	PULSE_STIM The pulse stim command can be used to control delivery of a small number of pulses.	57
static short	RESET_CYCLE Reset Cycle command can be used to turn stim on then off	57
static short	RESTORE_CLINICIAN_SETTINGS The restore clinician settings.	58
static short	SEND_TRIGGER The send trigger.	59
static short	SET_ACTIVE_GROUP The set active group.	57
static short	SET_NEXUS_CONFIGURATION The set nexus configuration.	56
static short	START_REALTIME The start realtime.	56
static short	START_SENSING The start sensing.	56

static short	STOP_REALTIME The stop realtime.	56
static short	STOP_SENSING The stop sensing.	56
static short	THERAPY_OFF The therapy off.	57
static short	THERAPY_ON The therapy on.	57

Field Detail

GET_INS_INFO

```
public static final short GET_INS_INFO
```

The get ins info.

SET_NEXUS_CONFIGURATION

```
public static final short SET_NEXUS_CONFIGURATION
```

The set nexus configuration.

START_SENSING

```
public static final short START_SENSING
```

The start sensing.

STOP_SENSING

```
public static final short STOP_SENSING
```

The stop sensing.

START_REALTIME

```
public static final short START_REALTIME
```

The start realtime.

STOP_REALTIME

```
public static final short STOP_REALTIME
```

The stop realtime.

GET_STATUS

```
public static final short GET_STATUS
```

The get status.

GET_REALTIME_DATA

```
public static final short GET_REALTIME_DATA
```

The get realtime data.

CYCLE_STIM

```
public static final short CYCLE_STIM
```

The cycle stim.

RESET_CYCLE

```
public static final short RESET_CYCLE
```

Reset Cycle command can be used to turn stim on then off

PULSE_STIM

```
public static final short PULSE_STIM
```

The pulse stim command can be used to control delivery of a small number of pulses.

THERAPY_ON

```
public static final short THERAPY_ON
```

The therapy on.

THERAPY_OFF

```
public static final short THERAPY_OFF
```

The therapy off.

SET_ACTIVE_GROUP

```
public static final short SET_ACTIVE_GROUP
```

The set active group.

INC_AMP

```
public static final short INC_AMP
```

The inc amp.

DEC_AMP

```
public static final short DEC_AMP
```

The dec amp.

INC_PULSEWIDTH

```
public static final short INC_PULSEWIDTH
```

The inc pulsewidth.

DEC_PULSEWIDTH

```
public static final short DEC_PULSEWIDTH
```

The dec pulsewidth.

INC_FREQ

```
public static final short INC_FREQ
```

The inc freq.

DEC_FREQ

```
public static final short DEC_FREQ
```

The dec freq.

RESTORE_CLINICIAN_SETTINGS

```
public static final short RESTORE_CLINICIAN_SETTINGS
```

The restore clinician settings.

SEND_TRIGGER

```
public static final short SEND_TRIGGER
```

The send trigger.

Class CommandResponseCodes

[mdt.neuro.nexus.commands](#)

[java.lang.Object](#)

└ [mdt.neuro.nexus.commands.CommandResponseCodes](#)

```
public class CommandResponseCodes
extends Object
```

The Class CommandResponseCodes.

Field Summary		Page
static short	AMP_AT_LOWER_LIM The amp at lower lim.	63
static short	AMP_AT_UPPER_LIM The amp at upper lim.	63
static short	CHANGE_NOT_DONE The change not done.	64
static short	COMMAND_NOT_VALID_FOR_DEVICE The command is not valid for the device type (i.e.	66
static short	DEVICE_RESET The device reset.	64
static short	FREQ_AT_LOWER_LIM The freq at lower lim.	64
static short	FREQ_AT_UPPER_LIM The freq at upper lim.	64
static short	GROUP_INVALID The group invalid.	62
static short	GROUP_OUT_OF_RANGE The group out of range.	62
static short	INS_FILTERING_ON INS Filtering On (i.e.	66
static short	INS_POR The ins por.	66
static short	INVALID_AMP The invalid amp.	64
static short	INVALID_CMD The invalid cmd.	65
static short	INVALID_CMD_DATA The invalid cmd data.	65
static short	INVALID_PROG_IN_SET The invalid program in the active set.	62
static short	INVALID_REALTIME_RESPONSE The invalid realtime response.	66
static short	NO_ACTIVE_SET_DEFINED The no active set defined.	63
static short	NO_APP_BACKUP The no app backup.	63

static short	NO_TELEM_RESP The no telem resp.	65
static short	NOT_ACTIVE_GROUP The not active group.	65
static short	NOT_ALL_CMDS_EXECUTED The not all cmds executed.	66
static short	NOT_READY The not ready.	65
static short	NULL_GROUP_ACTIVE The null group active.	65
static short	NULL_PACKETS The null packets.	66
static short	PIC_BUSY The pic busy.	67
static short	PIC_IF_BUSY The pic if busy.	67
static short	PIC_INVALID_ADDR The pic invalid addr.	68
static short	PIC_INVALID_CMD The pic invalid cmd.	67
static short	PIC_INVALID_CMD_TYPE The pic invalid cmd type.	67
static short	PIC_INVALID_LEN The pic invalid len.	67
static short	PIC_INVALID_PARAM The pic invalid param.	67
static short	PIC_NO_RESP The pic no resp.	67
static short	PIC_READ_BACK_FAIL The pic read back fail.	67
static short	PIC_TEST_IN_PROG The pic test in prog.	68
static short	PROG_NOT_IN_GROUP The prog not in group.	63
static short	PROG_OUT_OF_RANGE The prog out of range.	62
static short	PROGRAM_INVALID The program invalid.	62
static short	PW_AT_LOWER_LIM The pw at lower lim.	63
static short	PW_AT_UPPER_LIM The pw at upper lim.	63
static short	REALTIME_NOT_ACTIVE The realtime not active.	65
static short	SENSE_CONFIG_ERROR The sense config error.	65
static short	SET_TRANSITION_IN_PROG The set transition in prog.	64

static short	SUCCESS The success.	62
static short	TELEM_ERROR The telem error.	66
static short	THERAPY_NOT_ON The therapy not on.	63
static short	THERAPY_OOR The therapy oor.	64

Field Detail

SUCCESS

```
public static final short SUCCESS
```

The success. (0)

PROG_OUT_OF_RANGE

```
public static final short PROG_OUT_OF_RANGE
```

The prog out of range. (10)

GROUP_OUT_OF_RANGE

```
public static final short GROUP_OUT_OF_RANGE
```

The group out of range. (11)

INVALID_PROG_IN_SET

```
public static final short INVALID_PROG_IN_SET
```

The invalid program in the active set. (12)

PROGRAM_INVALID

```
public static final short PROGRAM_INVALID
```

The program invalid. (13)

GROUP_INVALID

```
public static final short GROUP_INVALID
```

The group invalid. (14)

THERAPY_NOT_ON

public static final short THERAPY_NOT_ON

The therapy not on. (15)

NO_APP_BACKUP

public static final short NO_APP_BACKUP

The no app backup. (18)

PROG_NOT_IN_GROUP

public static final short PROG_NOT_IN_GROUP

The prog not in group. (20)

NO_ACTIVE_SET_DEFINED

public static final short NO_ACTIVE_SET_DEFINED

The no active set defined. (21)

AMP_AT_UPPER_LIM

public static final short AMP_AT_UPPER_LIM

The amp at upper lim. (22)

AMP_AT_LOWER_LIM

public static final short AMP_AT_LOWER_LIM

The amp at lower lim. (23)

PW_AT_UPPER_LIM

public static final short PW_AT_UPPER_LIM

The pw at upper lim. (24)

PW_AT_LOWER_LIM

public static final short PW_AT_LOWER_LIM

The pw at lower lim. ([25](#))

FREQ_AT_UPPER_LIM

```
public static final short FREQ_AT_UPPER_LIM
```

The freq at upper lim. ([26](#))

FREQ_AT_LOWER_LIM

```
public static final short FREQ_AT_LOWER_LIM
```

The freq at lower lim. ([27](#))

THERAPY_OOR

```
public static final short THERAPY_OOR
```

The therapy oor. ([29](#))

CHANGE_NOT_DONE

```
public static final short CHANGE_NOT_DONE
```

The change not done. ([30](#))

DEVICE_RESET

```
public static final short DEVICE_RESET
```

The device reset. ([32](#))

SET_TRANSITION_IN_PROG

```
public static final short SET_TRANSITION_IN_PROG
```

The set transition in prog. ([50](#))

INVALID_AMP

```
public static final short INVALID_AMP
```

The invalid amp. ([51](#))

NOT_ACTIVE_GROUP

```
public static final short NOT_ACTIVE_GROUP
```

The not active group. (54)

NULL_GROUP_ACTIVE

```
public static final short NULL_GROUP_ACTIVE
```

The null group active. (59)

NOT_READY

```
public static final short NOT_READY
```

The not ready. (101)

INVALID_CMD

```
public static final short INVALID_CMD
```

The invalid cmd. (102)

INVALID_CMD_DATA

```
public static final short INVALID_CMD_DATA
```

The invalid cmd data. (103)

REALTIME_NOT_ACTIVE

```
public static final short REALTIME_NOT_ACTIVE
```

The realtime not active. (105)

SENSE_CONFIG_ERROR

```
public static final short SENSE_CONFIG_ERROR
```

The sense config error. (106)

NO_TELEM_RESP

```
public static final short NO_TELEM_RESP
```

The no telem resp. ([107](#))

TELEM_ERROR

```
public static final short TELEM_ERROR
```

The telem error. ([108](#))

INS_FILTERING_ON

```
public static final short INS_FILTERING_ON
```

INS Filtering On (i.e. real-time will not start if filtering is enabled in the Activa PC+S) ([109](#))

INS_POR

```
public static final short INS_POR
```

The ins por. ([110](#))

INVALID_REALTIME_RESPONSE

```
public static final short INVALID_REALTIME_RESPONSE
```

The invalid realtime response. ([111](#))

NULL_PACKETS

```
public static final short NULL_PACKETS
```

The null packets. ([112](#))

NOT_ALL_CMDS_EXECUTED

```
public static final short NOT_ALL_CMDS_EXECUTED
```

The not all cmds executed. ([114](#))

COMMAND_NOT_VALID_FOR_DEVICE

```
public static final short COMMAND_NOT_VALID_FOR_DEVICE
```

The command is not valid for the device type (i.e. sensing command for PC device). ([115](#))

PIC_IF_BUSY

```
public static final short PIC_IF_BUSY
```

The pic if busy. ([232](#))

PIC_INVALID_LEN

```
public static final short PIC_INVALID_LEN
```

The pic invalid len. ([233](#))

PIC_NO_RESP

```
public static final short PIC_NO_RESP
```

The pic no resp. ([234](#))

PIC_READ_BACK_FAIL

```
public static final short PIC_READ_BACK_FAIL
```

The pic read back fail. ([235](#))

PIC_BUSY

```
public static final short PIC_BUSY
```

The pic busy. ([241](#))

PIC_INVALID_CMD_TYPE

```
public static final short PIC_INVALID_CMD_TYPE
```

The pic invalid cmd type. ([242](#))

PIC_INVALID_PARAM

```
public static final short PIC_INVALID_PARAM
```

The pic invalid param. ([243](#))

PIC_INVALID_CMD

```
public static final short PIC_INVALID_CMD
```

The pic invalid cmd. ([244](#))

PIC_INVALID_ADDR

```
public static final short PIC_INVALID_ADDR
```

The pic invalid addr. ([245](#))

PIC_TEST_IN_PROG

```
public static final short PIC_TEST_IN_PROG
```

The pic test in prog. ([246](#))

Class CycleStim

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.CycleStim
```

```
public class CycleStim
extends ApplicationCommand
```

The Class CycleStim.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	69

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT

The timeout value (ms) for this command.
```

Class DataStructureAccessor

[mdt.neuro.nexus.commands](#)

[java.lang.Object](#)

└ [mdt.neuro.nexus.commands.DataStructureAccessor](#)

```
abstract public class DataStructureAccessor
extends Object
```

The Class DataStructureAccessor. This class allows friend access to the creation of data structures by internal API classes not in the data package.

Constructor Summary	Page
DataStructureAccessor()	70

Method Summary	Page
protected abstract DataPacket newDataPacket (List<Byte> bytes) Provides a New data packet.	71
protected abstract InsInfo newInsInfo (List<Byte> bytes) Provides a New ins info.	70
protected abstract NexusStatus newNexusStatus (List<Byte> bytes) Provides a New nexus status.	71
static void setProvider (DataStructureAccessor provider) Sets the provider.	70

Constructor Detail

DataStructureAccessor

```
public DataStructureAccessor()
```

Method Detail

setProvider

```
public static void setProvider(DataStructureAccessor provider)
```

Sets the provider.

Parameters:

provider - the new provider

newInsInfo

```
protected abstract InsInfo newInsInfo(List<Byte> bytes)
```

Provides a New ins info.

Parameters:

bytes - the bytes

Returns:

the ins info

newNexusStatus

```
protected abstract NexusStatus newNexusStatus(List<Byte> bytes)
```

Provides a New nexus status.

Parameters:

bytes - the bytes

Returns:

the nexus status

newDataPacket

```
protected abstract DataPacket newDataPacket(List<Byte> bytes)
```

Provides a New data packet.

Parameters:

bytes - the bytes

Returns:

the data packet

Class DecAmplitude

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├─ mdt.neuro.nexus.ApplicationCommand
│   └─ mdt.neuro.nexus.commands.DecAmplitude
```

```
public class DecAmplitude
extends ApplicationCommand
```

The Class DecAmplitude.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	72
static int	TIMEOUT	72

Fields inherited from class mdt.neuro.nexus.ApplicationCommand

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	73
double	getNewAmplitude () Gets the new amplitude.	73
byte	getNumExecuted () Gets the number of DEC commands executed – the actual number executed may not match the number requested.	73

Methods inherited from class mdt.neuro.nexus.ApplicationCommand

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

Method Detail

getNewAmplitude

```
public double getNewAmplitude()
```

Gets the new amplitude.

Returns:
the newAmplitude

getNumExecuted

```
public byte getNumExecuted()
```

Gets the number of DEC commands executed – the actual number executed may not match the number requested.

Returns:
the number of DEC commands executed

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class DecFrequency

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.DecFrequency
```

```
public class DecFrequency
extends ApplicationCommand
```

The Class DecFrequency.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	74
static int	TIMEOUT The timeout value (ms) for this command.	74

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes , commandCode , crcValid , request , response , responseCode , responseData , timeout

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	75
int	getNewFrequency () Gets the new frequency.	75

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
getRequestData , getResponseBytes , getResponseData

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getNewFrequency

```
public int getNewFrequency()
```

Gets the new frequency.

Returns:
the newFrequency

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class DecPulseWidth

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.DecPulseWidth
```

```
public class DecPulseWidth
    extends ApplicationCommand
```

The Class DecPulseWidth.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	76
static int	TIMEOUT The timeout value (ms) for this command.	76

Fields inherited from class mdt.neuro.nexus.ApplicationCommand

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	77
short	getNewPulseWidth () Gets the new pulse width.	77

Methods inherited from class mdt.neuro.nexus.ApplicationCommand

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getNewPulseWidth

```
public short getNewPulseWidth()
```

Gets the new pulse width.

Returns:
the newPulseWidth

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class GetInsInfo

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├─ mdt.neuro.nexus.ApplicationCommand
│   └─ mdt.neuro.nexus.commands.GetInsInfo
```

```
public class GetInsInfo
extends ApplicationCommand
```

The Class GetInsInfo.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	78
static int	TIMEOUT The timeout value (ms) for this command.	78

Fields inherited from class mdt.neuro.nexus.ApplicationCommand

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	79
InsInfo	getInfo () Gets the info.	79

Methods inherited from class mdt.neuro.nexus.ApplicationCommand

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getInfo

```
public InsInfo getInfo()
```

Gets the info.

Returns:
the info

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class GetNexusStatus

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.GetNexusStatus
```

```
public class GetNexusStatus
    extends ApplicationCommand
```

The Class GetNexusStatus.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	80

Fields inherited from class mdt.neuro.nexus. ApplicationCommand
commandBytes , commandCode , crcValid , request , response , responseCode , responseData , timeout

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	81
NexusStatus	getStatus () Gets the status.	80

Methods inherited from class mdt.neuro.nexus. ApplicationCommand
getRequestData , getResponseBytes , getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getStatus

```
public NexusStatus getStatus ()
```

Gets the status.

Returns:
the status

consumeBytes

protected void **consumeBytes** ([List<Byte>](#) bytes)

Overrides:

[consumeBytes](#) in class [ApplicationCommand](#)

Class GetRealtimeData

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.GetRealtimeData
```

```
public class GetRealtimeData
    extends ApplicationCommand
```

The Class GetRealtimeData.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	82

Fields inherited from class [mdt.neuro.nexus.ApplicationCommand](#)

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	83
DataPacket	getPacket () Gets the packet.	82

Methods inherited from class [mdt.neuro.nexus.ApplicationCommand](#)

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getPacket

```
public DataPacket getPacket ()
```

Gets the packet.

Returns:
the packet

consumeBytes

protected void **consumeBytes** ([List<Byte>](#) bytes)

Overrides:

[consumeBytes](#) in class [ApplicationCommand](#)

Class IncAmplitude

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.IncAmplitude
```

```
public class IncAmplitude
extends ApplicationCommand
```

The Class IncAmplitude.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	84
static int	TIMEOUT The timeout value (ms) for this command.	84

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes , commandCode , crcValid , request , response , responseCode , responseData , timeout

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	85
double	getNewAmplitude () Gets the new amplitude.	85
byte	getNumExecuted () Gets the number of INC commands executed – the actual number executed may not match the number requested.	85

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
getRequestData , getResponseBytes , getResponseData

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getNewAmplitude

```
public double getNewAmplitude()
```

Gets the new amplitude.

Returns:
the newAmplitude

getNumExecuted

```
public byte getNumExecuted()
```

Gets the number of INC commands executed – the actual number executed may not match the number requested.

Returns:
the number of INC commands executed

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class IncFrequency

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.IncFrequency
```

```
public class IncFrequency
extends ApplicationCommand
```

The Class IncFrequency.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	86
static int	TIMEOUT The timeout value (ms) for this command.	86

Fields inherited from class mdt.neuro.nexus.ApplicationCommand

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	87
int	getNewFrequency () Gets the new frequency.	87

Methods inherited from class mdt.neuro.nexus.ApplicationCommand

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getNewFrequency

```
public int getNewFrequency()
```

Gets the new frequency.

Returns:
the newFrequency

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class IncPulseWidth

mdt.neuro.nexus.commands

```

java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.IncPulseWidth

```

```

public class IncPulseWidth
extends ApplicationCommand

```

The Class IncPulseWidth.

Field Summary		Page
static int	RESPONSE_SIZE The response size.	88
static int	TIMEOUT The timeout value (ms) for this command.	88

Fields inherited from class mdt.neuro.nexus.ApplicationCommand

commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Method Summary		Page
protected void	consumeBytes (List<Byte> bytes)	89
short	getNewPulseWidth () Gets the new pulse width.	89

Methods inherited from class mdt.neuro.nexus.ApplicationCommand

getRequestData, getResponseBytes, getResponseData

Field Detail

RESPONSE_SIZE

```
public static final int RESPONSE_SIZE
```

The response size.

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Method Detail

getNewPulseWidth

```
public short getNewPulseWidth()
```

Gets the new pulse width.

Returns:
the newPulseWidth

consumeBytes

```
protected void consumeBytes(List<Byte> bytes)
```

Overrides:
[consumeBytes](#) in class [ApplicationCommand](#)

Class **PulseStim**

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├── mdt.neuro.nexus.ApplicationCommand
│   └── mdt.neuro.nexus.commands.PulseStim
```

```
public class PulseStim
    extends ApplicationCommand
```

The Class **PulseStim**.

Field Summary		Page
static int	responseSize The response size.	90
static int	timeout	90

Fields inherited from class [mdt.neuro.nexus.ApplicationCommand](#)

[commandBytes](#), [commandCode](#), [crcValid](#), [request](#), [response](#), [responseCode](#), [responseData](#), [timeout](#)

Constructor Summary		Page
PulseStim (byte progNumber, byte numSteps, byte numNop)	Instantiates a new pulse stim command.	91

Method Summary		Page
void	consumeBytes (List<Byte> bytes)	91
double	getNewAmplitude () Gets the new amplitude.	91
byte	getNumExecuted () Gets the number of DEC commands executed – the actual number executed may not match the number requested.	91

Methods inherited from class [mdt.neuro.nexus.ApplicationCommand](#)

[getRequestData](#), [getResponseBytes](#), [getResponseData](#)

Field Detail

responseSize

```
public static final int responseSize
```

The response size.

timeout

```
public static final int timeout
```

Constructor Detail

PulseStim

```
public PulseStim(byte progNumber,  
                 byte numSteps,  
                 byte numNop)
```

Instantiates a new pulse stim command.

Parameters:

progNumber - the prog number
numSteps - the num steps to INC then DEC
numNop - the num nop's to execute in between INC and DEC

Method Detail

getNewAmplitude

```
public double getNewAmplitude()
```

Gets the new amplitude.

Returns:

the newAmplitude

getNumExecuted

```
public byte getNumExecuted()
```

Gets the number of DEC commands executed – the actual number executed may not match the number requested.

Returns:

the number of DEC commands executed

consumeBytes


```
public void consumeBytes(List<Byte> bytes)
```

Overrides:

[consumeBytes](#) in class [ApplicationCommand](#)

Class `ResetCycle`

`mdt.neuro.nexus.commands`

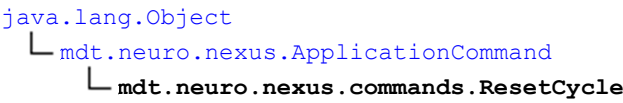

```
public class ResetCycle
extends ApplicationCommand
```

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	92

Fields inherited from class <code>mdt.neuro.nexus.ApplicationCommand</code>
<code>commandBytes</code> , <code>commandCode</code> , <code>crcValid</code> , <code>request</code> , <code>response</code> , <code>responseCode</code> , <code>responseData</code> , <code>timeout</code>

Methods inherited from class <code>mdt.neuro.nexus.ApplicationCommand</code>
<code>consumeBytes</code> , <code>getRequestData</code> , <code>getResponseBytes</code> , <code>getResponseData</code>

Field Detail

TIMEOUT

```
public static final int TIMEOUT

    The timeout value (ms) for this command.
```

Class RestoreClinicianSettings

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.RestoreClinicianSettings
```

```
public class RestoreClinicianSettings
extends ApplicationCommand
```

The Class RestoreClinicianSettings.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	93

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class SendTrigger

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.SendTrigger
```

```
public class SendTrigger
extends ApplicationCommand
```

The Class SendTrigger.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	94

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT

The timeout value (ms) for this command.
```

Class SetActiveGroup

mdt.neuro.nexus.commands

```
java.lang.Object
├─ mdt.neuro.nexus.ApplicationCommand
│   └─ mdt.neuro.nexus.commands.SetActiveGroup
```

```
public class SetActiveGroup
extends ApplicationCommand
```

The Class SetActiveGroup.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	95

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class SetNexusConfiguration

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.SetNexusConfiguration
```

```
public class SetNexusConfiguration
extends ApplicationCommand
```

The Class SetNexusConfiguration.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	96

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class StartRealTime

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.StartRealTime
```

```
public class StartRealTime
extends ApplicationCommand
```

The Class StartRealTime.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	97

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT

The timeout value (ms) for this command.
```

Class StartSensing

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.StartSensing
```

```
public class StartSensing
extends ApplicationCommand
```

The Class StartSensing.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	98

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class StopRealTime

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.StopRealTime
```

```
public class StopRealTime
extends ApplicationCommand
```

The Class StopRealTime.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	99

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class StopSensing

mdt.neuro.nexus.commands

```
java.lang.Object
└─ mdt.neuro.nexus.ApplicationCommand
    └─ mdt.neuro.nexus.commands.StopSensing
```

```
public class StopSensing
extends ApplicationCommand
```

The Class StopSensing.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	100

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class TherapyOff

[mdt.neuro.nexus.commands](#)

```
java.lang.Object
├─ mdt.neuro.nexus.ApplicationCommand
│   └─ mdt.neuro.nexus.commands.TherapyOff
```

```
public class TherapyOff
    extends ApplicationCommand
```

The Class TherapyOff.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	101

Fields inherited from class mdt.neuro.nexus. ApplicationCommand
commandBytes , commandCode , crcValid , request , response , responseCode , responseData , timeout

Methods inherited from class mdt.neuro.nexus. ApplicationCommand
consumeBytes , getRequestData , getResponseBytes , getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Class TherapyOn

mdt.neuro.nexus.commands

```
java.lang.Object
├─ mdt.neuro.nexus.ApplicationCommand
│   └─ mdt.neuro.nexus.commands.TherapyOn
```

```
public class TherapyOn
extends ApplicationCommand
```

The Class TherapyOn.

Field Summary		Page
static int	TIMEOUT The timeout value (ms) for this command.	102

Fields inherited from class mdt.neuro.nexus.ApplicationCommand
commandBytes, commandCode, crcValid, request, response, responseCode, responseData, timeout

Methods inherited from class mdt.neuro.nexus.ApplicationCommand
consumeBytes, getRequestData, getResponseBytes, getResponseData

Field Detail

TIMEOUT

```
public static final int TIMEOUT
```

The timeout value (ms) for this command.

Package mdt.neuro.nexus.data

Class Summary		Page
DataPacket	The Class DataPacket.	108
InsInfo	The Class InsInfo.	112
NexusStatus	The Class NexusStatus.	117
ProgramInfo	The Class ProgramInfo.	120

Enum Summary		Page
AmplitudeResolution	The Enum AmplitudeResolution.	103
ChannelType	The Enum ChannelType.	106
NexusState	The Enum NexusState.	115
SensingState	The Enum SensingState.	123
TherapyParameter	The Enum TherapyParameter.	125

Enum AmplitudeResolution

[mdt.neuro.nexus.data](#)

```
java.lang.Object
└─ java.lang.Enum<AmplitudeResolution>
    └─ mdt.neuro.nexus.data.AmplitudeResolution
```

All Implemented Interfaces:

[Comparable<AmplitudeResolution>](#), [Serializable](#)

```
public enum AmplitudeResolution
extends Enum<AmplitudeResolution>
```

The Enum AmplitudeResolution.

Enum Constant Summary	Page
CURRENT Constant current.	104
VOLTAGE Constant voltage.	104

Method Summary	Page
<small>static AmplitudeResolution</small> valueOf (String name)	105
<small>static AmplitudeResolution[]</small> values ()	104

Enum Constant Detail

VOLTAGE

```
public static final AmplitudeResolution VOLTAGE
```

Constant voltage.

CURRENT

```
public static final AmplitudeResolution CURRENT
```

Constant current.

Method Detail

values

```
public static AmplitudeResolution[] values()
```

valueOf

```
public static AmplitudeResolution valueOf(String name)
```

Enum ChannelType

[mdt.neuro.nexus.data](#)

```
java.lang.Object
└─ java.lang.Enum<ChannelType>
    └─ mdt.neuro.nexus.data.ChannelType
```

All Implemented Interfaces:

[Comparable<ChannelType>](#), [Serializable](#)

```
public enum ChannelType
extends Enum<ChannelType>
```

The Enum ChannelType.

Enum Constant Summary	Page
DISABLED The channel is disabled.	106
POWER The channel is power.	106
TIME The channel is time domain.	106

Method Summary	Page
<small>static ChannelType</small> valueOf (String name)	107
<small>static ChannelType[]</small> values ()	107

Enum Constant Detail

DISABLED

```
public static final ChannelType DISABLED
```

The channel is disabled.

POWER

```
public static final ChannelType POWER
```

The channel is power.

TIME

```
public static final ChannelType TIME
```

The channel is time domain.

Method Detail

values

```
public static ChannelType[] values()
```

valueOf

```
public static ChannelType valueOf(String name)
```

Class DataPacket

`mdt.neuro.nexus.data`

`java.lang.Object`

└ `mdt.neuro.nexus.data.DataPacket`

```
public class DataPacket
extends Object
```

The Class DataPacket.

Method Summary		Page
boolean	<code>convertToMillivolts(String fileName)</code> Converts the time domain data in this packet to mV, based on the given XML header file from the Sensing Programmer.	111
boolean	<code>equals(Object obj)</code>	111
byte	<code>getActiveGroup()</code> Gets the active group.	110
int[]	<code>getChSampleRates()</code> Gets the ch sample rates.	110
ChannelType[]	<code>getChTypes()</code> Gets the ch types.	110
double[][]	<code>getConvertedData()</code> Gets the converted data.	109
short[][]	<code>getData()</code> Gets the data.	110
short	<code>getNumMissedPatterns()</code> Gets the number of missed packets previous to this packet.	109
short	<code>getPatternNum1()</code> Gets the first pattern num.	109
short	<code>getPatternNum2()</code> Gets the second pattern num.	109
double[][]	<code>getSubtractedData()</code> Gets the template subtracted data.	108
int	<code>hashCode()</code>	111
byte	<code>isStimOn()</code> Checks if is stim on.	109
void	<code>subtractTemplate(double[][] template)</code> Subtract the given template from this packet and store the result in subtracted data.	110
String	<code>toString()</code>	111

Method Detail

getSubtractedData

```
public double[][] getSubtractedData()
```

Gets the template subtracted data.

Returns:
the template subtracted data

getConvertedData

```
public double[][] getConvertedData()
```

Gets the converted data.

Returns:
the converted data

getNumMissedPatterns

```
public short getNumMissedPatterns()
```

Gets the number of missed packets previous to this packet.

Returns:
the number of missed packets

getPatternNum1

```
public short getPatternNum1()
```

Gets the first pattern num.

Returns:
the first pattern num

getPatternNum2

```
public short getPatternNum2()
```

Gets the second pattern num.

Returns:
the second pattern num

isStimOn

```
public byte isStimOn()
```

Checks if is stim on.

Returns:
the stimOn

getActiveGroup

```
public byte getActiveGroup()
```

Gets the active group.

Returns:
the activeGroup

getData

```
public short[][] getData()
```

Gets the data.

Returns:
the data

getChSampleRates

```
public int[] getChSampleRates()
```

Gets the ch sample rates.

Returns:
the ch sample rates

getChTypes

```
public ChannelType[] getChTypes()
```

Gets the ch types.

Returns:
the ch types

subtractTemplate

```
public void subtractTemplate(double[][] template)
```

Subtract the given template from this packet and store the result in subtracted data.

Parameters:
template - the template

convertToMillivolts

```
public boolean convertToMillivolts(String fileName)
    throws ParserConfigurationException,
           org.xml.sax.SAXException,
           IOException
```

Converts the time domain data in this packet to mV, based on the given XML header file from the Sensing Programmer.

Parameters:

`fileName` - the file name of the XML header

Throws:

[ParserConfigurationException](#) - the parser configuration exception

[org.xml.sax.SAXException](#) - the SAX exception

[IOException](#) - Signals that an I/O exception has occurred.

hashCode

```
public int hashCode()
```

Overrides:

[hashCode](#) in class [Object](#)

equals

```
public boolean equals(Object obj)
```

Overrides:

[equals](#) in class [Object](#)

toString

```
public String toString()
```

Overrides:

[toString](#) in class [Object](#)

Class InsInfo

[mdt.neuro.nexus.data](#)

[java.lang.Object](#)

└─ [mdt.neuro.nexus.data.InsInfo](#)

```
public class InsInfo
    extends Object
```

The Class InsInfo. This class represents stimulation information and the status of sensing within the INS.

Field Summary		Page
static int	EXPECTED_SIZE The expected size.	112

Method Summary		Page
boolean	equals (Object obj)	113
int	getActiveGroupFrequency () Gets the active group frequency.	113
byte	getActiveGroupNumber () Gets the active group number.	113
ProgramInfo []	getPrograms () Gets the programs.	113
SensingState	getSensingState () Gets the sensing state.	113
int	hashCode ()	113
boolean	isTherapyOn () Checks if is therapy on.	112
String	toString ()	114

Field Detail

EXPECTED_SIZE

```
public static final int EXPECTED_SIZE
```

The expected size.

Method Detail

isTherapyOn

```
public boolean isTherapyOn()
```

Checks if is therapy on.

Returns:

the therapyOn

getSensingState

```
public SensingState getSensingState()
```

Gets the sensing state.

Returns:
the sensing state

getActiveGroupNumber

```
public byte getActiveGroupNumber()
```

Gets the active group number.

Returns:
the activeGroupNumber

getActiveGroupFrequency

```
public int getActiveGroupFrequency()
```

Gets the active group frequency.

Returns:
the activeGroupFrequency

getPrograms

```
public ProgramInfo[] getPrograms()
```

Gets the programs.

Returns:
the programs

hashCode

```
public int hashCode()
```

Overrides:
[hashCode](#) in class [Object](#)

equals

```
public boolean equals(Object obj)
```

Overrides:

`equals` in class `Object`

toString

`public String toString()`

Overrides:

`toString` in class `Object`

Enum NexusState

[mdt.neuro.nexus.data](#)

```
java.lang.Object
└─ java.lang.Enum<NexusState>
    └─ mdt.neuro.nexus.data.NexusState
```

All Implemented Interfaces:

[Comparable<NexusState>](#), [Serializable](#)

```
public enum NexusState
extends Enum<NexusState>
```

The Enum NexusState.

Enum Constant Summary	Page
IDLE The nexus is idle.	115
INS_CONNECTED The ins connected.	116
LINK_FAILED_DEVICE_ERR The link failed - device error.	116
LINK_FAILED_NO_RESPONSE The link failed - no response from INS.	116
LINKING_TO_INS The nexus is linking to INS.	115
MAINTENANCE_ENABLED The maintenance enabled.	116

Method Summary	Page
static NexusState valueOf (String name)	116
static NexusState[] values ()	116

Enum Constant Detail

IDLE

```
public static final NexusState IDLE
```

The nexus is idle.

LINKING_TO_INS

```
public static final NexusState LINKING_TO_INS
```

The nexus is linking to INS.

LINK_FAILED_NO_RESPONSE

```
public static final NexusState LINK_FAILED_NO_RESPONSE
```

The link failed - no response from INS.

LINK_FAILED_DEVICE_ERR

```
public static final NexusState LINK_FAILED_DEVICE_ERR
```

The link failed - device error.

INS_CONNECTED

```
public static final NexusState INS_CONNECTED
```

The ins connected.

MAINTENANCE_ENABLED

```
public static final NexusState MAINTENANCE_ENABLED
```

The maintenance enabled.

Method Detail

values

```
public static NexusState[] values()
```

valueOf

```
public static NexusState valueOf(String name)
```

Class NexusStatus

mdt.neuro.nexus.data

java.lang.Object

└ mdt.neuro.nexus.data.NexusStatus

```
public class NexusStatus
extends Object
```

The Class NexusStatus.

Field Summary		Page
static int	EXPECTED_SIZE The Constant EXPECTED_SIZE.	117

Method Summary		Page
boolean	equals (Object obj)	119
double	getBatteryPercent () Gets the battery percent.	118
byte	getHostTimeoutMinutes () Gets the host timeout minutes.	118
byte	getMaintenanceTimeoutSeconds () Gets the maintenance timeout seconds.	119
byte	getMajorVersion () Gets the major version.	118
byte	getMinorVersion () Gets the minor version.	118
NexusState	getState () Gets the state.	117
int	hashCode ()	119
boolean	isBatteryDepleted () Checks if is battery depleted.	118
String	toString ()	119

Field Detail

EXPECTED_SIZE

```
public static final int EXPECTED_SIZE
```

The Constant EXPECTED_SIZE.

Method Detail

getState

```
public NexusState getState()
```

Gets the state.

Returns:
the state

getMajorVersion

```
public byte getMajorVersion()
```

Gets the major version.

Returns:
the majorVersion

getMinorVersion

```
public byte getMinorVersion()
```

Gets the minor version.

Returns:
the minorVersion

getBatteryPercent

```
public double getBatteryPercent()
```

Gets the battery percent.

Returns:
the batteryPercent

isBatteryDepleted

```
public boolean isBatteryDepleted()
```

Checks if is battery depleted.

Returns:
the batteryDepleted

getHostTimeoutMinutes

```
public byte getHostTimeoutMinutes()
```

Gets the host timeout minutes.

Returns:
the hostTimeoutMinutes

getMaintenanceTimeoutSeconds

```
public byte getMaintenanceTimeoutSeconds()
```

Gets the maintenance timeout seconds.

Returns:
the maintenanceTimeoutSeconds

hashCode

```
public int hashCode()
```

Overrides:
[hashCode](#) in class [Object](#)

equals

```
public boolean equals(Object obj)
```

Overrides:
[equals](#) in class [Object](#)

toString

```
public String toString()
```

Overrides:
[toString](#) in class [Object](#)

Class ProgramInfo

[mdt.neuro.nexus.data](#)[java.lang.Object](#)└ [mdt.neuro.nexus.data.ProgramInfo](#)

```
public class ProgramInfo
extends Object
```

The Class ProgramInfo. This class represents the state of a stimulation program in the INS.

Field Summary		Page
static double	AMP_STEP_SIZE_CC The constant current amp step size in mA.	120
static double	AMP_STEP_SIZE_CV The constant voltage amp step size in V.	120

Method Summary		Page
boolean	equals (Object obj)	121
double	getAmplitude () Gets the amplitude in V.	121
AmplitudeResolution	getAmplitudeResolution () Gets the amplitude resolution.	121
byte	getProgramIndex () Gets the program index.	121
int	getPulseWidth () Gets the pulse width in uSec.	121
int	hashCode ()	121
String	toString ()	122

Field Detail

AMP_STEP_SIZE_CV

```
public static final double AMP_STEP_SIZE_CV
```

The constant voltage amp step size in V.

AMP_STEP_SIZE_CC

```
public static final double AMP_STEP_SIZE_CC
```

The constant current amp step size in mA.

Method Detail

getProgramIndex

```
public byte getProgramIndex()
```

Gets the program index.

Returns:
the programIndex

getAmplitude

```
public double getAmplitude()
```

Gets the amplitude in V. This amplitude is the actual value.

Returns:
the amplitude

getAmplitudeResolution

```
public AmplitudeResolution getAmplitudeResolution()
```

Gets the amplitude resolution. This value is set by the 8840 to be 1, which corresponds to 50 mV steps.

Returns:
the amplitudeResolution

getPulseWidth

```
public int getPulseWidth()
```

Gets the pulse width in uSec.

Returns:
the pulseWidth

hashCode

```
public int hashCode()
```

Overrides:
[hashCode](#) in class [Object](#)

equals

```
public boolean equals(Object obj)
```

Overrides:

`equals` in class `Object`

toString

`public String toString()`

Overrides:

`toString` in class `Object`

Enum SensingState

[mdt.neuro.nexus.data](#)

```
java.lang.Object
└─ java.lang.Enum<SensingState>
    └─ mdt.neuro.nexus.data.SensingState
```

All Implemented Interfaces:

[Comparable<SensingState>](#), [Serializable](#)

```
public enum SensingState
extends Enum<SensingState>
```

The Enum SensingState.

Enum Constant Summary	Page
DISABLED Sensing disabled.	123
ENABLED Sensing enabled.	123
UNAVAILABLE Status unavailable.	123

Method Summary	Page
<small>static SensingState</small> valueOf (String name)	124
<small>static SensingState[]</small> values ()	124

Enum Constant Detail

DISABLED

```
public static final SensingState DISABLED
```

Sensing disabled.

ENABLED

```
public static final SensingState ENABLED
```

Sensing enabled.

UNAVAILABLE

```
public static final SensingState UNAVAILABLE
```

Status unavailable.

Method Detail

values

```
public static SensingState[] values()
```

valueOf

```
public static SensingState valueOf(String name)
```

Enum TherapyParameter

[mdt.neuro.nexus.data](#)

```
java.lang.Object
└─ java.lang.Enum<TherapyParameter>
    └─ mdt.neuro.nexus.data.TherapyParameter
```

All Implemented Interfaces:

[Comparable<TherapyParameter>](#), [Serializable](#)

```
public enum TherapyParameter
extends Enum<TherapyParameter>
```

The Enum TherapyParameter.

Enum Constant Summary	Page
AMPLITUDE The Amplitude.	125
FREQUENCY The Frequency.	125
PULSEWIDTH The Pulse width.	125

Method Summary	Page
static TherapyParameter valueOf (String name)	126
static TherapyParameter[] values ()	126

Enum Constant Detail

AMPLITUDE

```
public static final TherapyParameter AMPLITUDE
```

The Amplitude.

PULSEWIDTH

```
public static final TherapyParameter PULSEWIDTH
```

The Pulse width.

FREQUENCY

```
public static final TherapyParameter FREQUENCY
```

The Frequency.

Method Detail

values

```
public static TherapyParameter[] values()
```

valueOf

```
public static TherapyParameter valueOf(String name)
```

Package mdt.neuro.nexus.support

Class Summary		Page
NexusLogger	The Class NexusLogger.	127

Class NexusLogger

[mdt.neuro.nexus.support](#)

`java.lang.Object`
└─ `mdt.neuro.nexus.support.NexusLogger`

```
public class NexusLogger
extends Object
```

The Class NexusLogger.

Constructor Summary	Page
NexusLogger()	128

Method Summary	Page		
<table><tr><td><code>static Logger</code></td><td>getLogger() Gets the logger.</td></tr></table>	<code>static Logger</code>	getLogger() Gets the logger.	128
<code>static Logger</code>	getLogger() Gets the logger.		

Constructor Detail

NexusLogger

```
public NexusLogger()
```

Method Detail

getLogger

```
public static Logger getLogger()
```

Gets the logger.

Returns:
the logger

